

Data Sheet

DIBROMO-TYROSINE (DiBrY)

ANTIBODY, MONOCLONAL

| | |
|----------------------------|--|
| Catalog no.: | AA1014.1 |
| Immunogen: | 3, 5-dibromo-tyrosine KLH-conjugated |
| Host: | Mouse Balb/c |
| Clone no.: | 3A5 |
| Isotype: | IgG ₁ kappa |
| Matrix: | Protein A purified, 10 mM PBS pH 7.4 |
| Specificity: | <p>Specific for di-halogenated tyrosine. Recognizes dibromo-tyrosine (DiBrY), N-acetylated dichloro-tyrosine (NAcDiClY), diiodo-tyrosine (DiIY), N-bromo-succinimide-treated BSA, HOCl-treated BSA and dibromo-tyrosine conjugated BSA. Bromotyrosine (BrY) is weakly recognized.</p> <p>There was no cross reactivity obtained with dityrosine-conjugated BSA, peroxyntirite-treated BSA, tyrosine (Y), chlorotyrosine(ClY), 3-3'-dityrosine (DiY), 3-aminotyrosine (NH₂Y), 3-nitrotyrosine (NO₂Y), phenylalanine (Phe), o-phosphotyrosine, 3, 4-dihydroxyphenylacetic acid (DOPA) and 2-hydroxy-4-phenylacetic acid. No cross reactivity was obtained with brominated nucleic acids 5-bromouracil, 5-bromouridine, and 8-bromoguanosine.</p> |
| Contents: | <p>20 µg (lyophilized)</p> <p>Resuspend in 200 µl aqua bidest. for 100 µg/ml of IgG</p> |
| Known applications: | ELISA ¹ , Western Blot, immunohistochemistry (paraffin sections, 10 µg/ml) ¹ |



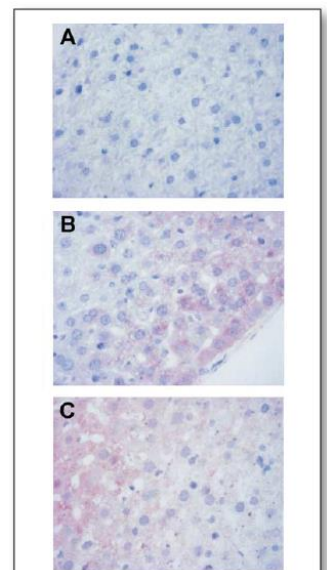


Figure 1: Immunohistochemistry image of Dibromo-tyrosine staining in paraffin sections of LPS-treated mouse liver. Antigen retrieval was performed by microwave treatment in 0.01 M citrate buffer (pH 6.0) for 5 min. The sections were incubated with AA1014 (10 µg/ml) and detected using biotinylated secondary antibody and Vectastain ABC-AP kit (Vector Laboratories). Vector Red substrate served as chromogen (Vector Laboratories). Sections were counterstained with hematoxylin. **A.** Liver section of untreated mouse; **B.** liver section of LPS-treated mouse sacrificed after 24 h; **C.** liver section of LPS-treated mouse sacrificed after 72 h.

Kato Y et al. (2005) *Free Radical Biology & Medicine* 38:24–31

This antibody has not been tested for use in all applications. This does not necessarily exclude its use in non-tested procedures. The stated dilutions are recommendations only. End users should determine optimal dilutions in their system using appropriate negative/positive controls.

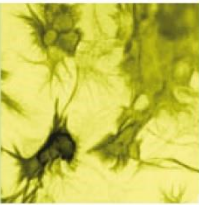

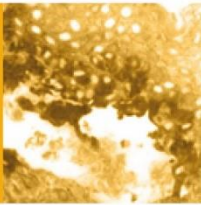

Store at: - 20 °C

Repeated thawing and freezing must be avoided

References:

1. Kato Y, Kawai Y, Morinaga H, Kondo H, Dozaki N, Kitamoto N, Osawa T. (2005). Immunogenicity of a brominated protein and successive establishment of a monoclonal antibody to dihalogenated tyrosine. *Free Radic Biol Med* **38**(1): 24-31
2. Wu W, Chen Y, d'Avignon A, Hazen SL. (1999). 3-Bromotyrosine and 3,5-dibromotyrosine are major products of protein oxidation by eosinophil peroxidase: potential markers for eosinophil-dependent tissue injury in vivo. *Biochemistry* **38**(12):3538-48

Last updated on: 14 April 2022

For research use only

Publishing research using AA1014? Please let us know so that we can cite your publication as a reference.

