



Data Sheet


OXYTOCIN

ANTIBODY, POLYCLONAL

Catalog no.:	AS1005.1 / AS1005.2
Immunogen:	Synthetic Oxytocin, bovine thyroglobulin conjugated
Swiss-Prot No:	P01175
Gene Information:	Gene Name: OXT GeneID: 280888
Host:	Rabbit
Matrix:	Serum
Specificity:	Human, bovine and other species Oxytocin No cross reaction was obtained with Lys-/Arg-vasopressin, neurophysin, or pituitary hormones (e.g. prolactin)
Contents:	20 µl / 100 µl (lyophilized) Resuspend in 20 µl / 100 µl aqua bidest.
Known applications:	RIA (1: 80 000) ^{1, 2, 3, 4, 5} , Western Blot, immunohistochemistry (paraffin sections, 1:500) ⁵ This antibody has not been tested for use in all applications. This does not necessarily exclude its use for non-tested procedures. The stated dilutions are recommendations only. We suggest that the applicant titrates the antibody in his/her system using appropriate negative/positive controls.
Store at:	2-8 °C (lyophilized); - 20 °C (dissolved) Repeated thawing and freezing must be avoided
References:	1. Schams D, Schmidt-Polex B, Kruse V (1979). Oxytocin determination by radioimmunoassay in cattle. I. Method and preliminary physiological data. <i>Acta Endocrinol</i> 92 (2):258-270. 2. Schams D, Lahlou-Kassi A, Glatzel P (1982). Oxytocin concentrations in peripheral blood during the oestrous cycle and after ovariectomy in two breeds of sheep with low and high fecundity. <i>J Endocrinol</i> 92 (1):9-13. 3. Schams D (1983). Oxytocin determination by radioimmunoassay. III. Improvement to subpicogram sensitivity and application to blood levels in cyclic cattle. <i>Acta Endocrinol</i> 103 (2):180-183. 4. Walters DL, Schams D, Schallenberger E (1984). Pulsatile secretion of gonadotrophins, ovarian steroids and ovarian oxytocin during the luteal phase of the oestrous cycle in the cow. <i>J Reprod Fertil</i> 71 (2):479-491. 5. Kruij TA, Vullings HG, Schams D, Jonis J, Klarenbeek A (1985). Immunocytochemical demonstration of oxytocin in bovine ovarian tissues. <i>Acta Endocrinol</i> 109 (4):537-542.
Last updated on:	28 April 2022



Antibodies

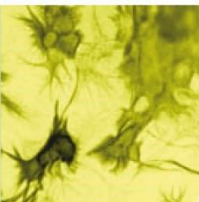

For research use only

Publishing research using AS1005? Please let us know so that we can cite your publication as a reference.



Immundiagnostik AG

Stubenwald-Allee 8a · 64625 Bensheim · Germany

Phone: +49 6251 70190-0 · Fax: +49 6251 70190-363 · dept.immuochemicals@immundiagnostik.com · www.immundiagnostik.com